

# UROMED STONIZER *advanced* Steinfangkörbchen

## UROMED STONIZER *advanced* Stone Retrieval Basket

**REF 76191212**

spitzenlos  
4-Draht-Körbchen, gerade


*tipless*  
**4-wire basket, straight**

### Der STONIZER *advanced*

Mit kreisrunder Korbform und Schlaufentechnik.

### The STONIZER *advanced*

*With circular basket design and loop technique.*

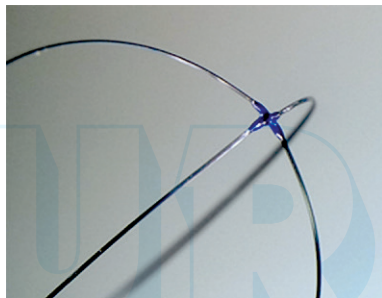

**Steinfangkörbchen zur urologischen Entfernung von Steinen aus dem Ureter und den Nieren während der Ureterorenoskopie**

- zur Anwendung in der Niere geeignet
- hohe Formstabilität durch Korbelement aus Nitinol
- Einhandbedienung
- Minikorb: im Vergleich zu REF 77191212 um 15% weniger raumfordernd (Korblänge 17,5mm)
- neuartige Konstruktion des Korbes und modifizierte Materialien sichern höchste Zuverlässigkeit – auch beim Bergen von Steinstraßen
- Zug- und Druckstabilität des Körbchens durch beständiges Material
- Mechanismus zur Begrenzung der Zugkraft gewährleistet sichere Steinentfernung
- atraumatisches distales Korbeende
- Steuerung des Körbchens mit einer hohen Drehstabilität
- hohe Funktionalität auch bei abgewinkeltem Schaft
- flexibel durch Nitinol-Material
- aufgrund des geringen Außendurchmessers ist ein erhöhter Spülfluss möglich
- abnehmbarer Handgriff durch Rändelschraube
- zur einmaligen Verwendung
- Einführhilfe

**Retrieval basket for the urological removal of stones from the ureter and kidney during ureterorenoscopy**

- *for use in kidneys*
- *high shape stability due to nitinol basket*
- *one-handed use*
- *mini basket: 15% less space required compared to REF 77191212 (basket length 17,5mm)*
- *innovative construction of basket and modified materials ensure the highest possible reliability – also in case of stone streets*
- *pull and pressure stability of the basket because of its durable material*
- *mechanism to limit the tensile strength provides secure stone retrieval*
- *atraumatic distal end of basket*
- *guidance of the basket with high torque stability*
- *high functionality even with bended shaft*
- *flexible nitinol material*
- *increased flushing is possible due to the small outer diameter*
- *removable handle by knurled screw*
- *single use*
- *insertion guide*

REF	Ch./Fr.	Ø Schaft	Länge	Ø Korb
REF	Ch./Fr.	Ø Shaft	Length	Ø Basket
76191212	Ch. 1,9	0,6 mm	1200 mm	12 mm

**Verpackung:** steril  
**Originalkarton:** 1 Stück

**Packaging:** sterile  
**Original carton:** 1 piece